

# Mighty Light PLUS – Tungsten Halogen (MLP-TH)

*A powerful source of low-noise white light with versatile output modules*

The Mighty Light PLUS (MLP) is a broadband light source that provides 10X higher spatial brightness than competitive sources: delivering up to 7 W of collimated output from a 10 mm diameter flexible light guide. Useful output spans 300 to 2500 nm and the low-noise output power is smoothly adjustable from 0 - 100 %. Applications include microscopy, white light interferometry, machine vision, and precision inspection.

## Specifications


- Type of lamp : Tungsten Halogen
- Spectral wavelength : 300 - 2500 nm
- Lamp power consumption : 250 W
- Power stability : < 0.7 %
- Bulb lifetime : ~ 1000 hours
- Color temperature : 3,200 K
- Dimensions : 340 x 160 x 140 mm
- Power control knob : 0 - 100 %
- Electrical requirements : AC 100 - 240 V, 50/60 Hz




## Optional Specifications

- Correlated color temperature (CCT) : 3,400 K
- Bulb lifetime : ~ 500 hours

## Features

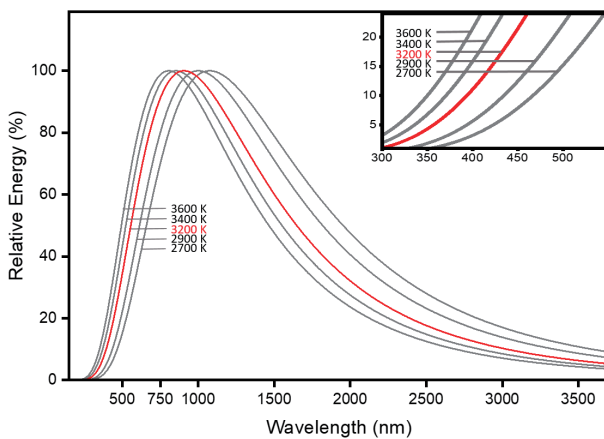

Broad output spectrum

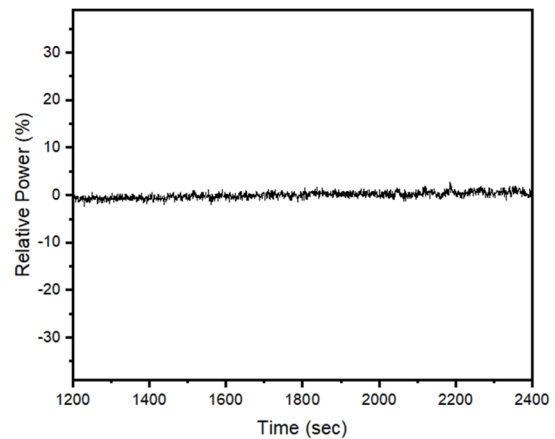

Low power fluctuation (High stability) : < 0.7 %